

## Reliable, confident, impressive.

A completely automatic management for an optimized procedure.  
With an outstanding simplicity.



Millseal Plus Automatic is an automatic unit that allows the operator to perform the whole procedure, starting from the bags preparation to their sealing, by means of simple keys.

	The motorized sliding in both directions allows to exactly define the bag length without paper wasting.
	A single key for a complete automatism of the sealing and cut steps.
	With a single touch, it is possible to seal up to 4 bags simultaneously.
	The program keys allows to store 2 different bag lengths and perform the sealing and cut steps automatically.



Technical Data	Millseal Plus Automatic	Millseal Plus Manual
Power supply	220/240 Vac	220/240 Vac
Mains frequency	50/60 Hz	50/60 Hz
Nominal power	150 W - 0.45 A	150 W - 0.45 A
External dimensions (without roll)	474 x 374 x 200 mm (LxDxH)	474 x 374 x 200 mm (LxDxH)
Weight	6 kg	5 kg
Type of sealing	With continuous cycle and constant temperature	With continuous cycle and constant temperature
Nominal operating temperature of the heating resistance	180 °C	180 °C
Height of the sealing band	12 mm	12 mm
Maximum sealing width	300 mm	300 mm
Sealable material	Paper/propylene film	Paper/propylene film
Maximum roll diameter	200 mm	200 mm

**mocom**

Via Saliceto, 22 - 40013 Castel Maggiore (BO) Italy - Tel. (+39) 051.700877 - Fax (+39) 051.701056  
mocomcom@mocom.it - [www.mocom.it](http://www.mocom.it)



Our thermosealers are manufactured in compliance with the following:

- Directives: 2006/95/CEE, 2004/108/CEE.
- Standards: EN 61010-1:2001, EN 61326-1:2006

Mocom S.r.l. reserves the right to make changes to any of the products shown in this document, without prior notice. VHAZ7080008



Millseal thermosealers



The top class thermosealers.

**mocom**

Naturally sterilizers

## Millseal Plus Automatic



## They protect and guarantee your professional competence.


High performance thermosealers, with an innovative design and a forefront technology. Indispensable.



### A guarantee for the sterilization maintenance

Millseal thermosealers are equipped with an electronic control of the sealing parameters (time and temperature) and assure a sealing band of 12 mm - so exceeding the requirements of the Norms.

The acoustic and visual signals provide information about the proper operative conditions of the thermosealing device.

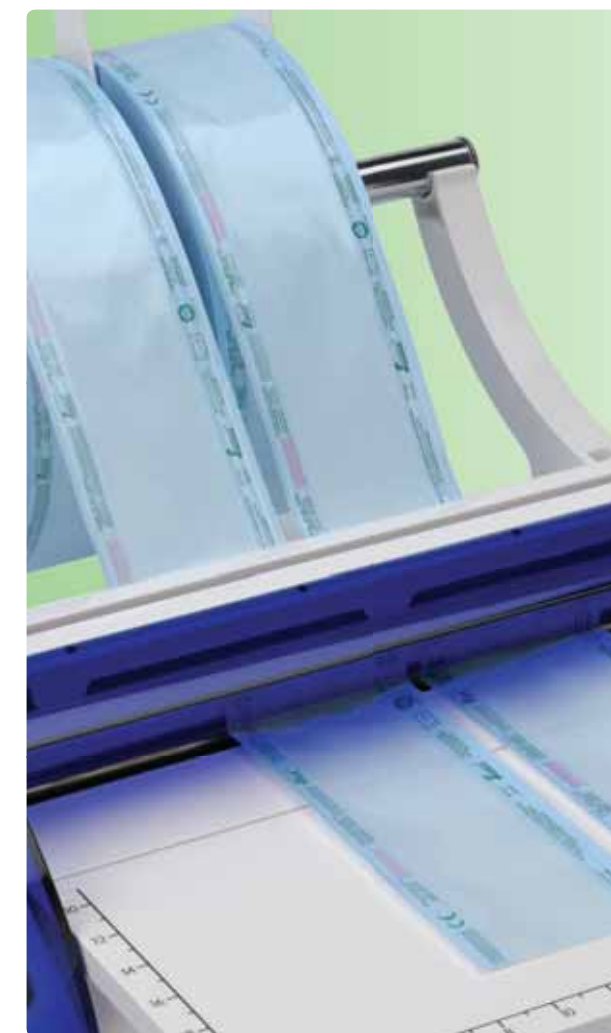

### High and steady performances

Thanks to the automatic adjustment of the temperature, it is possible to achieve an effective control of the sealing area heating, thus assuring high quality performances over time. Moreover, Millseal thermosealers are equipped with a safety device against overheating and with an automatic stand-by mode for energy saving.

### Ergonomic and functional design

Millseal thermosealers combine their ease of use with an outstanding reliability. The position maintenance system of the rolls avoid their re-winding after the cut and sealing steps, allowing a consistent and uninterrupted procedure.

The metric scale on the resting plate helps in defining with precision the proper length of the bags based on the type of material to be sterilized.

## Millseal Plus Manual

